



# Study Models Impression Checklist

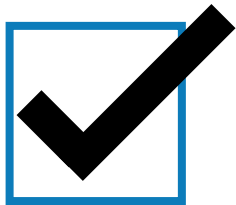

## Use a 5 day alginate

Almost every problem with bad impressions can be solved by using a 5 day alginate. For recommended alginates, please consult our 'Recommended Study Models Products List' document.\*

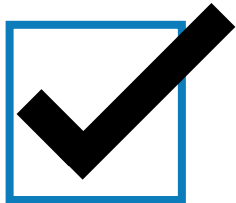

## Take a good impression

For tips on what constitutes a good impression, consult our 'Achieving Quality Impressions' document.\*



## Check impression

Once you take the impression, have a fellow team member or doctor double check the impression and make sure everything looks good.



## Take an impression bite

In addition to taking the impression, we ask that you also take a wax bite. A wax bite is important because it helps us get the teeth in accurate occlusion. For recommended wax bites, please consult our 'Recommended Study Models Products' document.\*




## Wrap impressions in a damp paper towel

Once impression is taken, wrap in a slightly damp paper towel to preserve the impression. Place covered impressions in a sealable plastic bag and pour a cap full of alcohol in the bag to prevent mold and bacteria from growing.



## Attach intraoral photos

In addition to the wax bite, intraoral photos also help us when determining if the study model matches the patient's teeth.

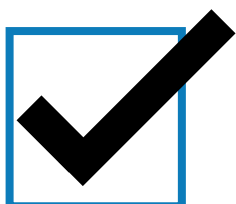

## Ship the case to OrthoSelect

Once all the previous steps are completed, place impressions in a box with sufficient padding and seal it. We pour up the impressions as soon as we receive them, so please ship the impressions within 2 days of when you take them. To qualify for free shipping, please include 2 sets of impressions in small box or 7 or more sets of impressions in a big box. For full instructions, please consult our 'Shipping Procedures' document.\*

\* - Found on our 'Study Models Resource Library' on [www.myorthoselect.com](http://www.myorthoselect.com) under the 'About' section